



MWB041

OPERATION CONNECTING TROLLEY

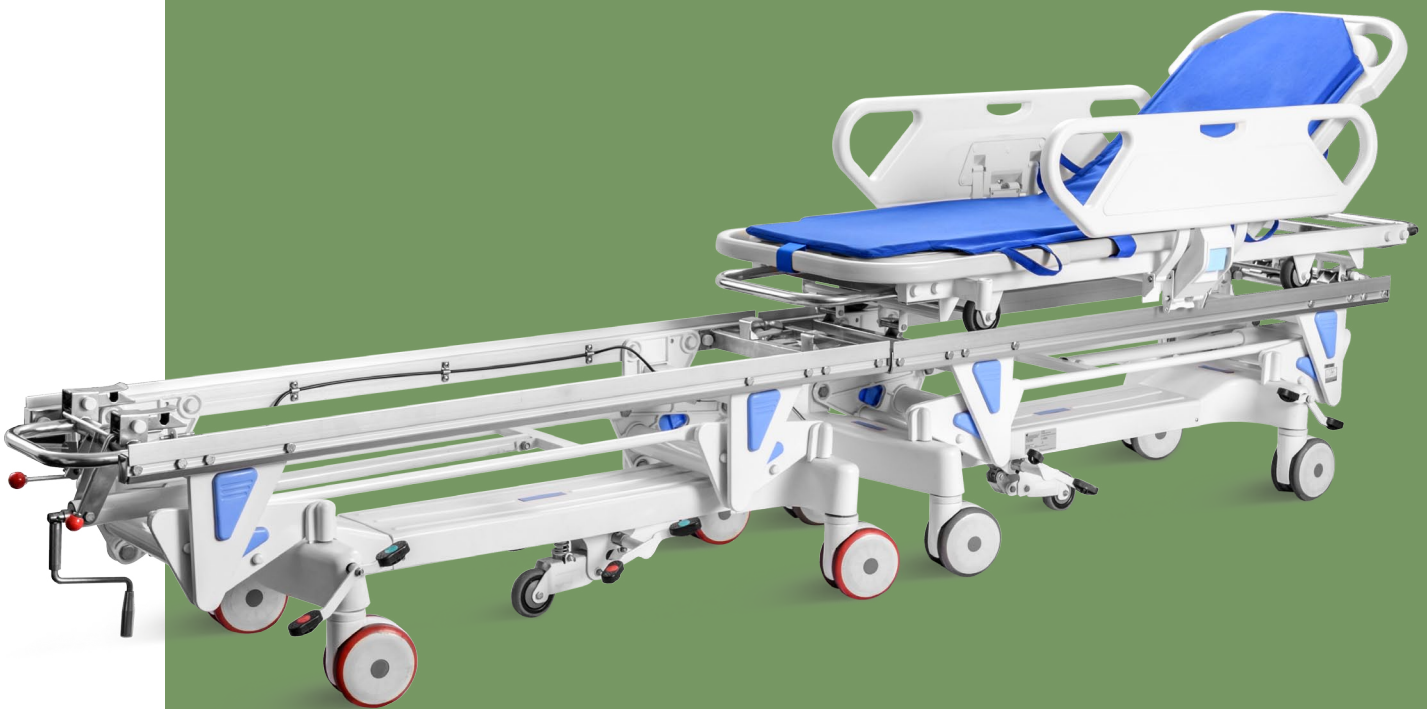

NWB041

FEATURES AND HIGHLIGHTS

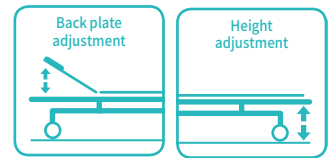


MAIN TECHNICAL PARAMETERS

▪ Size	3650*640*(640-970)mm
▪ Back-rest	0-75°(±10°)

TECHNICAL CONFIGURATION

☑ Crank System	2PC
☑ Back Rest Air Spring	IPC
☑ ABS Side Rail	1set
☑ Directional Wheel	2pc
☑ 6" Castors with Central Locking	4pcs
☑ IV Previsions	2pc
☑ IV Pole	1set
☑ Central Locking Break Pedal	4pcs
☑ Directional Wheel Release Pedal	4pcs

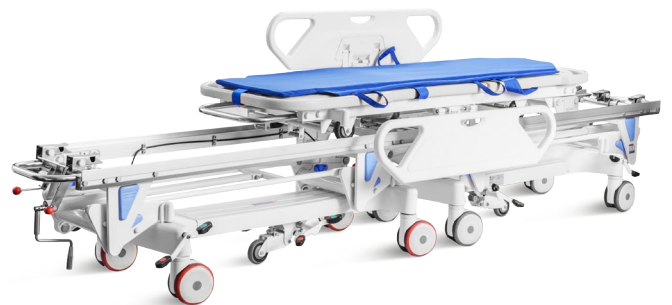



protect

cover all metal screws with rubber cap, protect screws from possible corrosion




Center Wheel

Can be engaged to facilitate direction control, enabling users to navigate the stretcher easily in hectic and dynamic environments.



connect with operating room

Transfer patient to operation room, prevent cross infection

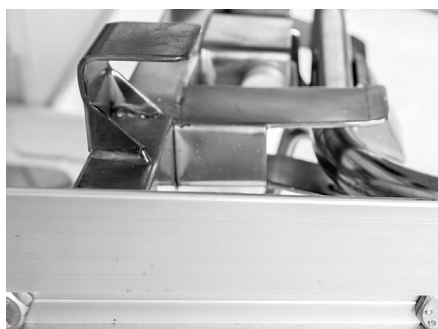

Metal Material

Baosteel from Fortune Global 500
Panasonic Robotic ensure 360° full smooth welding



Plastic

Environmental plastic can be sterilized in 100°C, tensile upto 30MPa



Steels | painting

11 Process epoxy painting, ASTM testing anti-bacterial, paint thickness 0.12mm, brightness 60°, paint can resist 50kg impact



Metal Dust Cover

We use metal dust protection for crank system, more durable than plastic cover



Bearing

We use robust tube in weight bearing structure



Autolock

The patient stretcher will be autolocked, when it slid to one side of the trolley



NUMBER WIN

MEDICAL TECHNOLOGY



7×24

Online 24 Hours for Sales Service
Timely Response for Clients' Problems

+86 18021231901

Suzhou Number Win Medical Technology Co.,Ltd.

No. 35, Lehong Road · Zhangjiagang (Jiangsu) · China

export@number-win.com

www.number-win.com

As the circumstance of printing, the photos might have minor difference on color, tactile from real products.
As the update of design and improvement, parameters and outlook of products might change, thanks for understanding.



Company video



General catalogue